


RifeMedic ® 7

www.rifemedic.co.za Rev 03 Dec.2020

Dual Rife Machine, Quantum Radionics Bio-Resonance Education and Application System

The RifeMedic 7 is a modern touch screen LCD multi-modality device with Educational features, Diagnosis support, Assessment and Scanner modalities, Bio-resonance applications, Health modality applications, and Tools for supporting a holistic approach to health care.

Two RifeMedic 7 LCD sizes are available: 7-inch ; 10-inch. Two RifeMedic 7 devices are offered: A double rife machine ; A double rife machine with scanner functions.

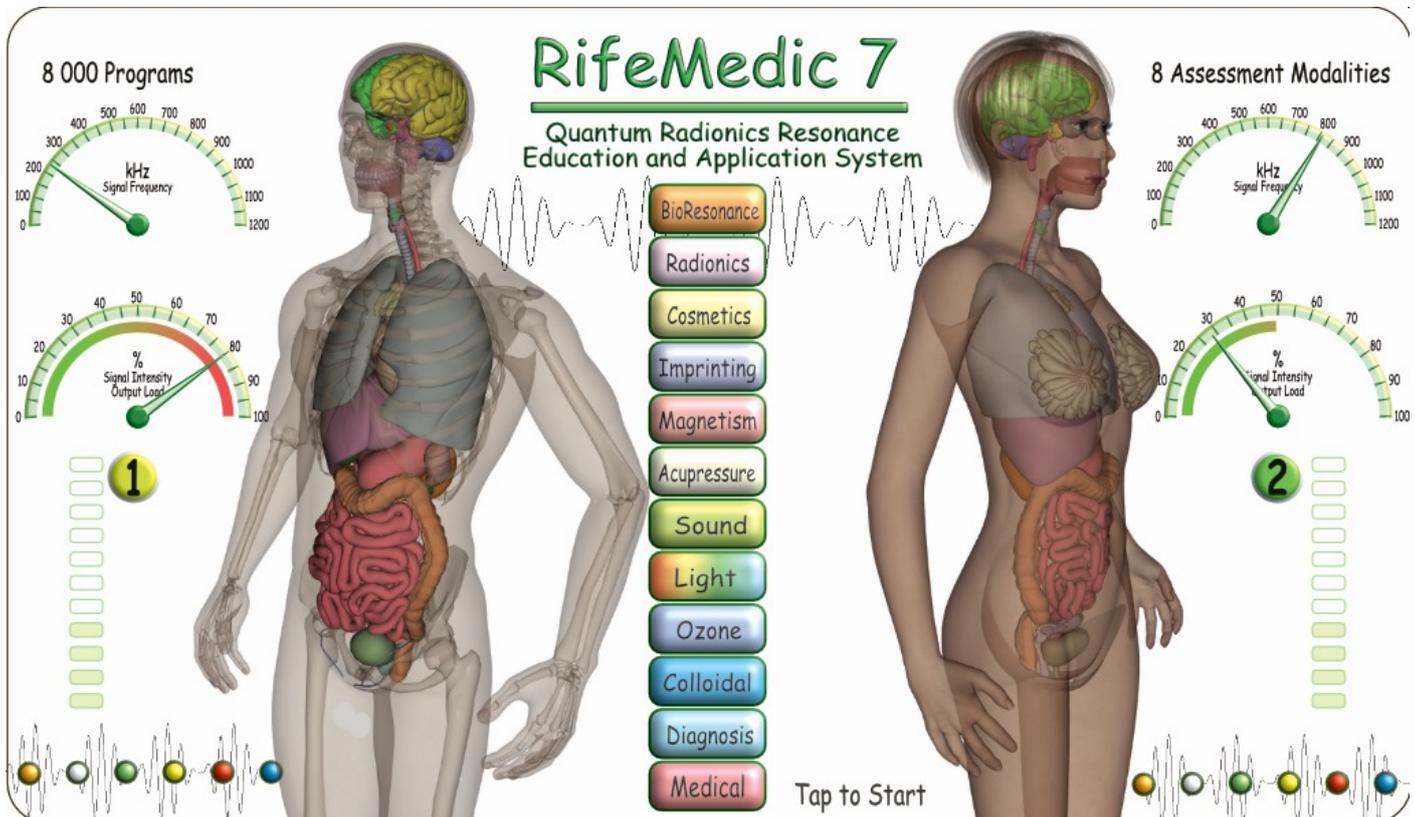

- RifeMedic 7 2XR-7 - 7" lcd dual bio-resonator
- RifeMedic 7 2XR-10 - 10" lcd dual bio-resonator
- RifeMedic 7 2XRA-7 - 7" lcd dual bio-resonator with assessment features and devices
- RifeMedic 7 2XRA-10 - 10" lcd dual bio-resonator with assessment features and devices

Features and functions of RifeMedic 7 devices are constantly being improved and updated.

Unique features include the Beat modes, Modulation modes and Position modes that have been unique features of the RifeMedic 5 devices. These modes apply dual frequency and triple frequency applications that have interestingly been recently co-discovered by other bio-resonance researchers to be extremely efficient and can not be compared with single frequency applications.

Specifications:

- 8 000 + pre-programmed rife machine programs
- 2 independent rife machine output channels
- 1 radionics and water imprint output channel (*) Note
- 8 health assessment and scanner modalities (*) Note
- 6 health support treatment modalities (*) Note
- 3 health related tools (*) Note
- 0 - 1 200 000 Hz bio-resonance frequency bandwidth per output channel
- 0 - 4 000 000 Hz bio-resonance frequency modulation system
- 48 V peak to peak maximum adjustable output per channel
- 15 nSec output switching speed

(*) Note: These features depend on the RifeMedic 7 model or may require extra devices to function.


Peripherals included: 2 x sets of hand electrodes, 1 x set of foot plates, 2 sets of connecting wires

Available peripherals: Electro-magnetic blanket, Radionics unit, Water-imprinter unit


Assessment peripherals and Tools: To be announced when ready.

Quick Start:

On all screens, press on the Man for a help screen

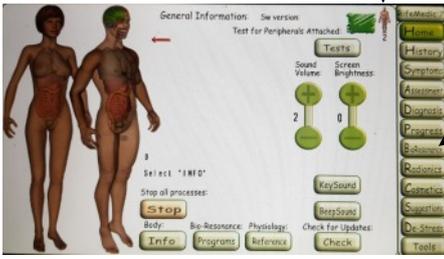

1: Press BioResonance
2: Press Programs

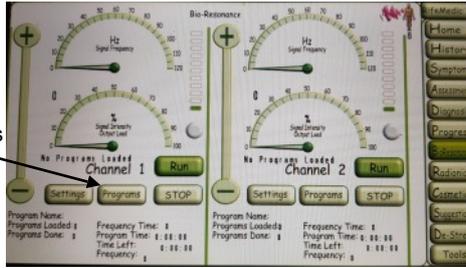


Programs Press Positions:

ON TOP is the CURRENTLY SELECTED PROGRAM

Increase the number
Decrease the number


(healing) Press the word to Search the next same word

- Go 1 Up in the list, direction = Up
- Go to Previous alpha letter, direction = Up
- Go 10 Up or Down the list
- Go to Next alpha letter, direction = Down
- Go 1 Down in the list, direction = Down

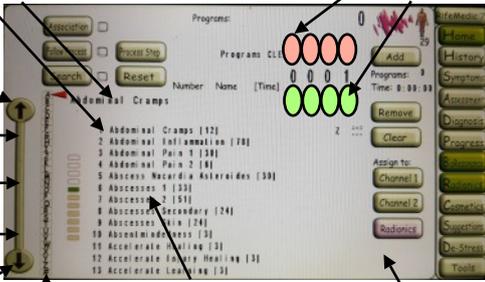

Select this Item
Select Segments of the Alphabet
Down one Page

Search

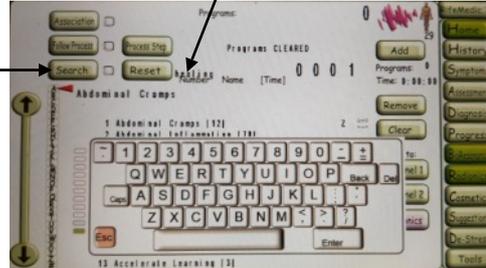


Type =healing to search for the word within the names
Select another list segment to search further from there
(go to next alpha letter and press the search word to search from there)

NB On the Programs Page: Do not make selections too fast to avoid screen sync loss - do not press wildly - take it calm

The top item is the currently selected program

3: Select and Assign Programs to one or both of the channels

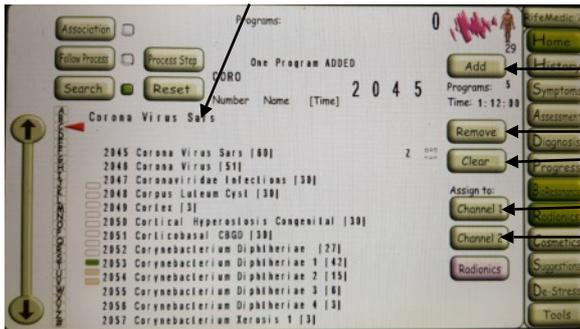


- Add the Program to the selection - the list of programs you wish to run
- Remove the last added program
- Clear the entire program selection
- Assign the selected programs to Channel 1
- Assign the selected programs to Channel 2

The procedure for assigning programs to Channel 1 and Channel 2 will be:
Add one or a selection of programs. Assign the loaded programs to Channel 1.
Clear the selection of programs.
Add one or a selection of programs. Assign the loaded programs to Channel 2.
Go to the BioResonance Page.

To delete programs after assigning it to a Channel, go to the Bio-Resonance Page And press the specific Channel's Stop button Twice.

Intensity select press points Automatic scale change

Remove the last added program

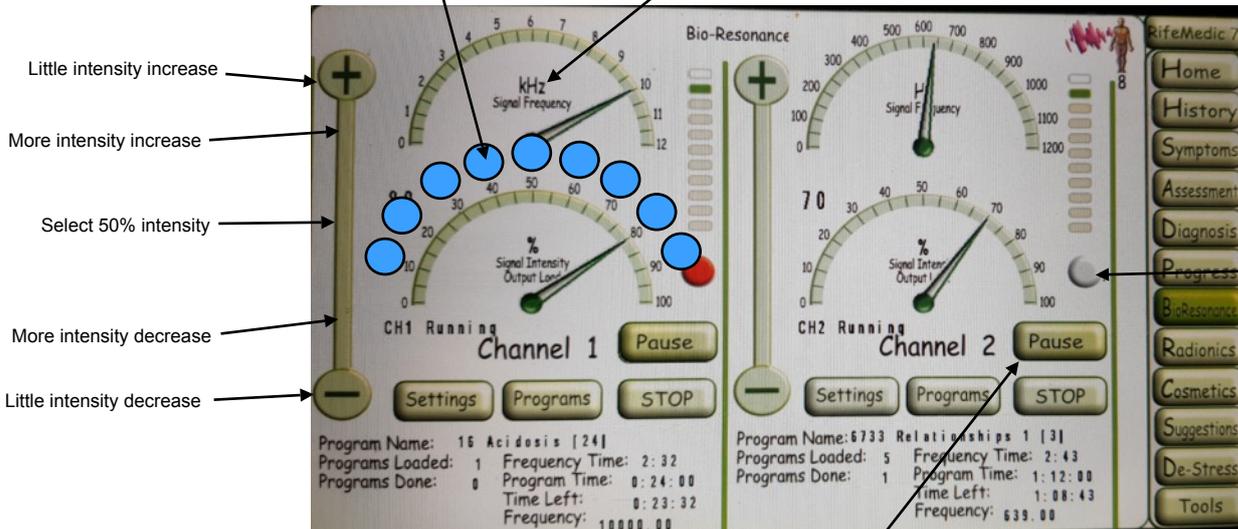

White = no programs loaded
Green = programs loaded
Pink = paused
Red = running

This button will be Run - press Run to run the Channel
Press Pause to pause the Channel - the Intensity selection will be reset

4: **While holding the electrodes in both hands or applied to both feet** - press Run and adjust the output to the desired intensity. Adjust the intensity until the sensation feels gentle and pleasant. Press intensity decrease any time while intensity increases to stop increasing or lower the intensity. Always check the output leds - they are directly connected to the electrode wires. High frequencies may not produce and sensation - do not adjust the output too high with high frequencies. Do not short out the electrodes.

5: In order to use the "Beat Mode" or the "Times Mode":

Beat Mode adds the selected frequency value of one channel to the frequency of the other channel.

For example, when you load the same programs in Channel 1 and Channel 2 like the exercise below, and you select "Beat Mode" on Channel 1 until 528 Hz appears, Channel 1's frequency = Channel 2's frequency + 528.


Times Mode multiplies the selected frequency value of one channel to the frequency of the other channel.

For example, when you load the same programs in Channel 1 and Channel 2 like the exercise below, and you select "Times Mode" on Channel 1 until 11 appears, Channel 1's frequency = Channel 2's frequency x 11.

Use only Beat Mode or Times Mode not both.

An example to let Channel 2 run 11 times faster than Channel 1:

1. Press "BioResonance" to go to the Bio-resonance page.
2. Press "Programs" to go to the Programs page
3. Press "Clear", then select your required programs for example "General Comprehensive" and press "Add" and then "Channel 1" and "Channel 2" so the same program is assigned to Channel 1 as well as to Channel 2.
4. Select the "BioResonance" page and press "Run" for Channel 1 and "Run" for Channel 2.
5. Select "Settings" and on channel 2, keep selecting "Times Mode" until the value is "11 x Frequency of Channel 1".
6. Select the "BioResonance" page and adjust the intensity of channel 1 to the point where it feels comfortable. Then adjust the intensity of channel 2 to the point where it feels comfortable. Observe the frequency of channel 2 is 11 times the frequency of channel 1.

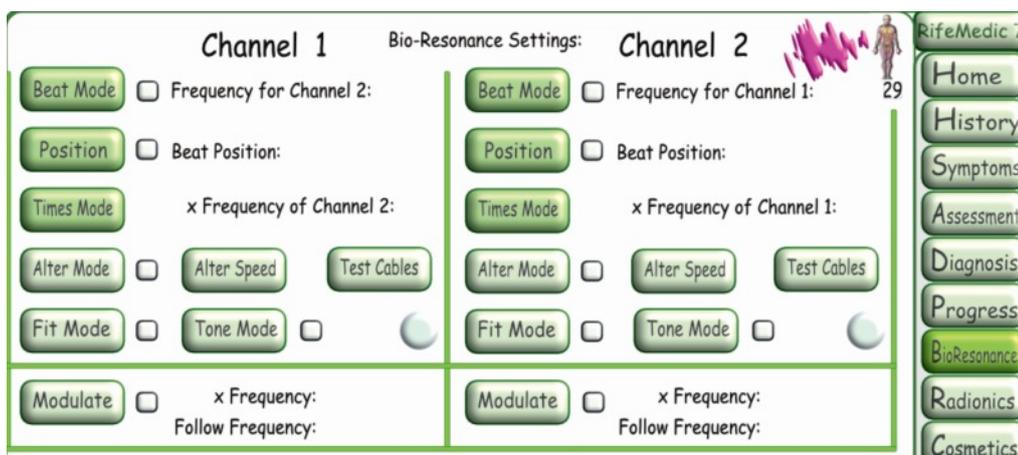

The RifeMedic 7 is a fully functional Bio-Resonance application device. It has many advanced features including a second channel that can function fully independent. Also different modes of operation and functioning for enhanced effect and other modality applications.

The RifeMedic 7 (2XR models) are development in-progress devices with new features being released about every week. Not all screens on these devices are active or usable as they may be for future implementation or only related to the more advanced scanner-RifeMedic 7 (2XRA models).

New in-progress features:

Alter Mode: Beat Mode and Times Mode Channel's frequencies are Altered or swapped

Modulation Modes: A modulation frequency can ride onto each Channel's output signal. Such a modulation frequency may enhance the currently applied frequency with an important high-frequency signal to further enhance the effect of the applied frequency.

Diagnosis Support: List your symptoms and the system will provide possible options in conditions that may associate with the symptoms.

Physiology Support: Several human physiology pages are available with vital information. Next will follow condition's effects and symptoms on the body.